

Powered by **Discover Next**

Move beyond the words in new ways

Incredible sound performance to hear the deeper meaning in conversations. Our latest platform: Discover Next. Now with an expanded product lineup and remote adjust.

Discover Next overview 2.0

A Sonova brand

unitron

It's all on Discover Next

Technology to hear beyond
sounds and words

Lithium-ion rechargeables
proven to last all day

A refreshed set of
easy-to-use accessories

Stylish hearing instruments
with all-day comfort

Easy connections to all
kinds of favorite tech

Enhanced flexibility
with remote fine-tuning
adjustments

Spectacular speech understanding is just the beginning

The sound technology in Discover Next delivers more than just sounds and words. It helps your clients hear the subtle nuances of speech by enabling better understanding of who is talking, where the voices are coming from, what intonation they're using, and what their underlying emotions are* – whether they're in a noisy environment or conversing with soft-spoken people.

We added more to love


An expanded product lineup featuring three new Inera™ ITE hearing instruments and Stride™ M BTE.

Increased flexibility in how you service clients with remote adjust – a new option that helps you deliver fine-tuning adjustments remotely.

*Singh, Kreuger, Besser, Wietoska, Launer, Meis. (2018). A Pre-Post Intervention Study of Hearing Aid Amplification: Results of the Emotional Communication in Hearing Questionnaire (EMO-CHeQ), IHCON 2018.

Help your clients get to the heart of conversations

SoundCore™ – the unique engine inside Discover Next – delivers the key nuances of speech and environmental cues that enables your clients to understand the deeper meaning in conversations.

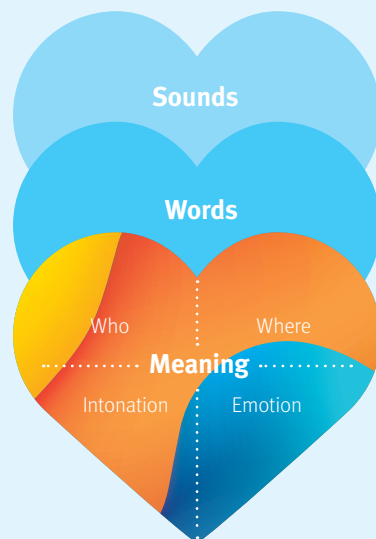


How we're moving beyond the words

Identifies environments as accurately as a 20-year-old

Clients can hear their best across all listening situations with **SoundNav 3.0** and can enjoy life-enriching media with the help of MediaNav. Years of training our algorithms with machine learning means this feature can identify life's ever-changing soundscapes just like a young, normal-hearing listener would.*

More voices. Less noise.

Provide clear speech, reduced noise, and overall sound quality that's comfortably balanced in any environment with **Sound Conductor**. And make the subtle nuances of speech more pronounced so everyone – even soft-spoken people – are easier to understand with our new **soft speech lift** feature.

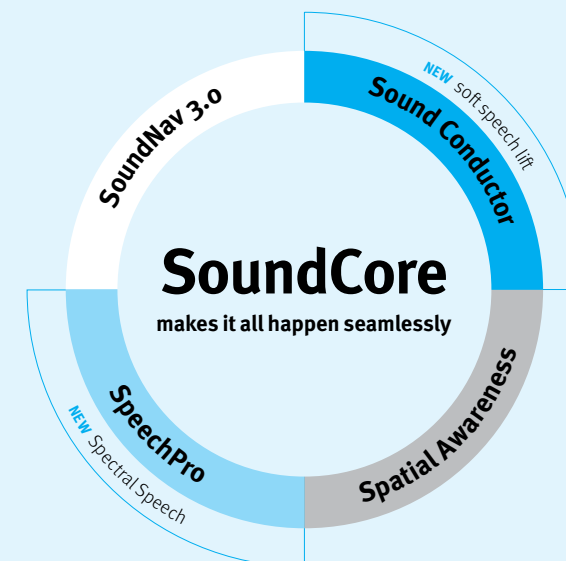

Hear speech with amazing clarity – even in noise

Address your clients' number one need to hear better in noisy environments with **SpeechPro**.** It uses binaural cues to know where speech is coming from and enhances it beyond what a typical hearing instrument can do. A new fourth dimension of SpeechPro, **Spectral Speech**, was designed to help clients better understand words and emotions in complex listening situations.

Put sounds in their place

Clients enjoy a realistic sound experience because sounds can be acoustically separated and perceived in their accurate locations. This is thanks to **Spatial Awareness**, which brings back natural localization cues with an innovative 4 microphone strategy.

*Dr. David A. Eddins, WCA October 2018 Cape Town, South Africa. **Ozmeral, Hoover, Gabbidon, and Eddins (2018). Development of Continuous Number Identification Test (CNIT). Ear and Hearing submitted for publication.

Proven to last and last and last

Say hello to our new family of stylish, comfortable hearing instruments, including an ITE with made for all (MFA) connectivity, and lithium-ion rechargeable options proven to last all day.



Moxi™ Move R
 Our smallest lithium-ion rechargeable. Receiver-in-canal.

Moxi Jump R T
 Rechargeable. Receiver-in-canal. Telecoil.

Moxi Fit
 Traditional battery. Receiver-in-canal.

NEW Stride M
 Our smallest behind-the-ear on Discover Next. Traditional battery. Telecoil.

Stride P R
 Rechargeable. Behind-the-ear.



NEW Insera W 312 Directional
 Our first in-the-ear with MFA connectivity. Traditional battery.

NEW Insera 312 Omni
 Traditional battery. In-the-ear. Telecoil.

NEW Insera 10A Omni
 Our smallest in-the-ear on Discover Next. Traditional battery. Telecoil.

Available colors

- Cinnamon (Q9)

Amber suede (S2)

Espresso (P4)

Platinum (P6)

Amber (P2)

Charcoal (P8)

Beige (01)

Pewter (P7)
- Available with ITEs
- Beige (26)

Cocoa (22)

Tan (14)

Brown (28)

Black* (06)

Hearing instruments shown at actual size.
 *Black is an extra option on the IIC style only.

Beyond world-class sound performance



Helpful solutions

Your clients will love the experience even more with our refreshed set of easy-to-use solutions.



Flexibility for you.

Convenience for all.

Life's easier for everyone with remote adjust

Move beyond the traditional way of supporting your clients with new **remote adjust**.

Enhanced flexibility is yours with the option to perform fine-tuning adjustments remotely. Remote adjust offers you the ability to deliver professional care that's convenient for you and your clients, whether they're at home or on the move.


Adjustments are easy to make and sent directly to your clients via the Remote Plus app*. With remote adjust, you can provide a hearing experience that's even easier and more empowering.



Clients can enjoy professional care without having to go into your clinic! In just a few taps on the Remote Plus app, they can apply your fine-tuning adjustments and hear the improved sound.

With the Remote Plus app your clients can also:


- › Receive helpful tips, advice and reminders on their smartphone that reinforce and extend your counseling with **coach**
- › Make adjustments to bass, mid and treble frequencies that are captured in our TrueFit™ software with the **equalizer**
- › Adjust the balance between media and the environment
- › Submit ratings on their experiences so you can be proactive with your counseling
- › Adjust volume, mute and unmute, and watch how-to videos

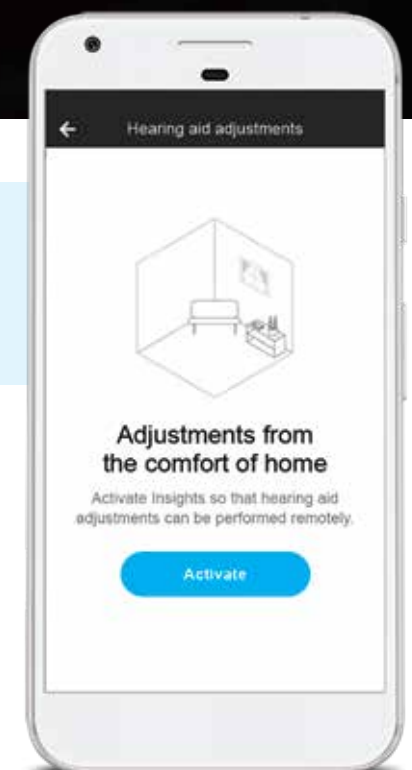


*For the use of the Remote Plus app, Unitron hearing aids with Bluetooth® connectivity are required. The Remote Plus app works on Apple smartphones with iOS 12 or newer and Android smartphones with version 7 or newer. App Store badge is a trademark of Apple Inc., registered in the U.S. and other countries. Google Play and the Google Play logo are trademarks of Google LLC.

Discover Next fitting ranges


Insera

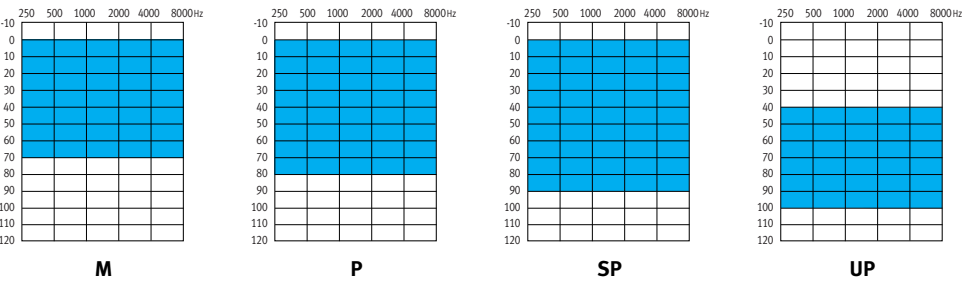

Stride

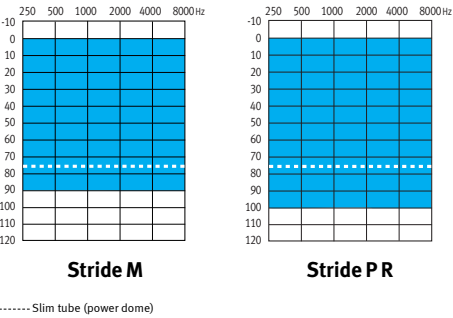

Moxi

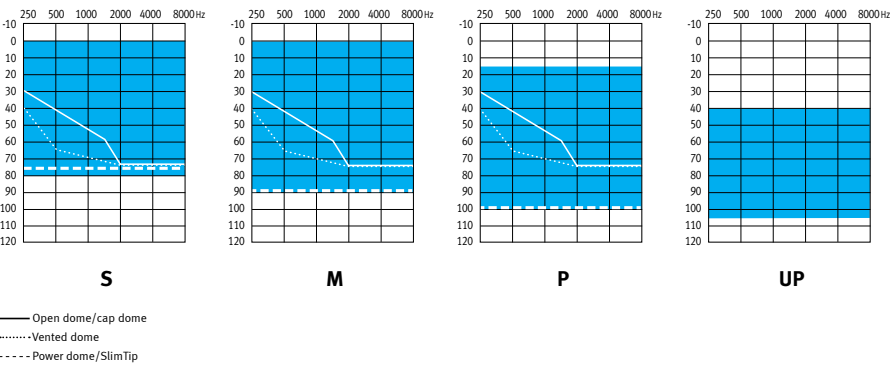

At a glance

| SoundCore | | Discover Next 9 | Discover Next 7 | Discover Next 5 | Discover Next 3 |
|---|-------------------------------|-----------------|-----------------|-----------------|-----------------|
| Environmental classification | SoundNav 3.0 | • | • | • | • |
| | Total listening environments | 7 | 6 | 4 | 2 |
| | Conversation in a crowd | • | | | |
| | Conversation in a small group | • | • | | |
| | Music | • | • | | |
| | Quiet | • | • | • | |
| | Noise | • | • | • | |
| | Conversation in quiet | • | • | • | • |
| | Conversation in noise | • | • | • | • |
| | Total streaming environments | 2 | 2 | 2 | 2 |
| Sound optimization | MediaNav speech | • | • | • | • |
| | MediaNav music | • | • | • | • |
| | Sound Conductor | • | • | • | • |
| | Speech enhancement | • | • | • | • |
| | Soft speech lift | • | • | • | |
| Localization | Noise reduction | • | • | • | • |
| | Microphone features | • | • | • | • |
| | Spatial Awareness | • | • | • | |
| | Personalized | • | | | |
| Performance in challenging environments | Pinna Effect | • | • | • | • |
| | SpeechPro | • | | | |
| | Speech Locator | • | • | | |
| | Speech Focus | • | | | |
| | Spatial Speech | • | | | |
| | Spectral Speech | • | | | |
| SpeechZone 2 | | | • | | |
| Fine-tuning channels | | 20 | 16 | 14 | 12 |

Available in all technology levels

| Sound stabilization | Ease and convenience | Experience innovations |
|---|---|---|
| <div><div>›</div> AntiShock 2</div> <div><div>›</div> Wind control</div> <div><div>›</div> Feedback manager</div> <div><div>›</div> Natural Sound Balance</div> | <div><div>›</div> Remote adjust</div> <div><div>›</div> Made for all direct connectivity</div> <div><div>›</div> TV Connector</div> <div><div>›</div> Remote Control</div> <div><div>›</div> PartnerMic</div> <div><div>›</div> Remote Plus app</div> <div><div>-</div> Equalizer</div> <div><div>›</div> DuoLink</div> <div><div>›</div> Binaural Phone</div> <div><div>›</div> RogerDirect™</div> <div><div>-</div> Roger Pen</div> <div><div>-</div> Roger Select</div> <div><div>-</div> Roger Table Mic II</div> | <div><div>›</div> FLEX:TRIAL™</div> <div><div>›</div> FLEX:UPGRADE™</div> <div><div>›</div> Insights</div> <div><div>-</div> Coach</div> <div><div>-</div> Log It All</div> <div><div>-</div> Ratings (including proactive alerts)</div> <div><div>-</div> Overall satisfaction</div> <div><div>-</div> Usage</div> |
| <div><div>›</div> Fitting</div> <div><div>›</div> First fit workflow</div> <div><div>›</div> Automatic Adaptation Manager</div> <div><div>›</div> Frequency compression 2</div> <div><div>›</div> Tinnitus masker</div> <div><div>›</div> MyMusic</div> <div><div>›</div> Manual programs</div> <div><div>›</div> IntelliVent</div> | | |

Compatibility may be limited by product style.

Even more to love



Easy and accurate first fits in 5 clicks or less, with sound that eases them in and keeps getting better.



Made for all connectivity for easy, direct connections to all the people and content they love.



An expanded product lineup based on award-winning designs, and the ability to perform professional fine-tuning adjustments remotely.

Differentiate yourself with FLEX™

From letting clients try hearing instruments before they buy, to upgrading technology in the same instrument and a whole lot more, FLEX is a game-changing approach that enhances the way you deliver hearing care.

To learn more, contact your representative today.



Unitron is a hearing solution company that believes people should feel really good about the entire hearing care experience, start to finish. Our ingenious products, technologies, services and programs offer a level of personalization you can't get anywhere else. Get ready to Love the experience.

© 2020 Unitron. All rights reserved.

A/2020-07 028-6459-02

unitron.com

sonova
HEAR THE WORLD